

Features:

Ionic detox footbath with dual system

Dual large LCD display

Dual advanced FIR beam

Strong steel main machine house

Aluminum + black leather carrying case

Benifits:

body waste & heavy metals removal

Oxygen & Energy Increased

Sleep improved

Weight loss

zon1.se

0736709495

2009 Latest Model

Luxurious Ionic Detox Footbath



*Ion Detox Therapy
User Manual*



• Disclaimer:

The system does not claim to cure or treat any disease or injury. It assists your body to re-balance its bio-energy field and stimulates the body for self-detoxification. The body's organs will naturally function better when the electro-magnetic fields are charged. The ion energy charged. Footbath will help you feel young again.

• The work theory of the ionic detox footbath

The ionic detox footbath process is an amazing cleansing method for detoxifying your entire body by bathing your feet in warm water while using positive and negative ions from an ion generator---the array. The process is also known as electrolysis. The array element is the actual ion generator device that sits down into the water and creates the ion charge and flow within that solution. It attaches to the power supply, which controls the array element ion flow and change.

• Benefits for Our Body:

- ◆ Removes toxins and body waste products that can cause health problems
- ◆ Helps clear up bad skin.
- ◆ Helps inactivate Viruses, Bacteria, Yeast, and Fungus
- ◆ Positive effect on weight loss (Assists in weight loss)
- ◆ Headache Relief - Assists in clearing up headache problems.
- ◆ Slows down aging and improves body flexibility.
- ◆ Rejuvenates and energizes your whole body.
- ◆ Enhances Nutrient Absorption
- ◆ Heavy Metals Removal
- ◆ Increase Your Energy
- ◆ Help with Faster Recovery Time for Disease and Injury
- ◆ Increased Oxygen in Your Body
- ◆ Helps Relieves Pain and Tension
- ◆ Improved Sleep
- ◆ Reduces Inflammation and Unwanted Fluid Retention
- ◆ Immune System Booster (Helps Balance Immune System)
- ◆ Safe for Children 8 and up

• Usage and Settings

The ionic detox footbath can be used by anyone over years of age, as long as the session is supervised and all instructions are followed

carefully. However, the recommended time and frequency is going to vary from person to person.

Children's sessions take considerably less than adults' sessions.

• General guidelines:

Adults under 50 years maximum 30 minute sessions, once every two days for 10-14 sessions

Adults over 50 years maximum 30 minute sessions, once every three days for 10-14 sessions

Children under 10 years 15 minute or per session, once every three days for 6-10 sessions

Use cautious judgment as everyone is different; ensure a two week break between programs.

Increase experiencing fatigue require angry, moody, or in a state of panic etc., allow more time between sessions and check for vitamin and mineral depletion. Fatigue indicates an overstressed condition, and emotional disturbance indicates that emotional toxin is being released along with chemical toxins.

• Colors and Particles in Water

Even if not bath the feet in water, the water also will change color because there is the dirt and chemical substances in water which can induce color change. Our body has 70% water which is absolutely necessary in life. The toxin produced from our body is very similar as from water. But when bathing feet, the dirt and toxin from water will not adhere to the container or edges of the electrode plate. The toxin from the body will form mucus supernatant that can be cleaned by the bleaching water or cleanser.

Many different colors and objects will appear in the water during the Ionic Detox Foot Bath sessions. Colors such as Yellow-green, orange, brown, black, dark green, white foam, white cheese like particles, black flecks, red flecks bubbles of differing colors, brown particles. Each color and object, as best as we have been able to determine based on correlation performed with an Ionic Detox Foot Bath.

The below can be for your reference after many tests:

- ◆ Yellow=green detoxifying from the kidney, bladder, urinary tract, female/prostate
- ◆ Orange= from joints
- ◆ Brown= liver, tobacco, cellular debris
- ◆ Black= liver, gallbladder

- ◆ Dark Green= gallbladder
- ◆ White foam= mucous from lymph
- ◆ White speckles= most likely yeast
- ◆ Black speckles= heavy metals
- ◆ Red speckles= blood clot material

The combination of the water, metal and salt will produce different colors and some substances even if not bath the feet in water. The theory of the Ionic Detox Foot Bath in water also is suitable for the body. But if different areas, the chemical pollution is different so with the different results. There are different popular symptoms in every area, such as bad arthritis, bladder. Suggestion is to test muscle after cleaning feet. The metal wastage of the electrode plate will create the metal substances in water. More and more experience and observes can help to identify where the toxin mainly is from. The stronger energy of the mineral in water will produce the more activate of Ionic Detox Foot Bath. If see from the result with many black spots, 20%-40% of which are impacted by local water quality, other 60%-80% are toxins from the body. If extremely producing bitter elements from the Ionic Detox Foot Bath machine, the water will not be suitable for use. If finding green color after steady test, it is the best to change distilled or mineral water with half spoon of sea salts or normal salts.

• Components of the system:

- ◆ Main machine with double LCD Display
- ◆ Array X 2 (10 coils) last 50-60 times
- ◆ New Double Infrared FIR Bamboo Charcoal Belt
- ◆ Power cable
- ◆ User Manual
- ◆ Dual wristbands
- ◆ Aluminum & black leather carrying case.
- ◆ 1 PH test kit
- ◆ 2 detox foot patches

• The ionic detox footbath unit (Main machine)



Black model




Silver model



Array



Case

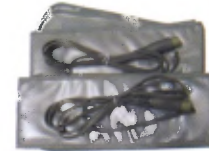

New Far Infrared belts



Power cord



Wristbands


• Features:

1. Large dual LCD display (9cm x 4.5cm). It is unparalleled in global market.
2. Show current curved line on the screen.
3. Solid carrying case (Aluminum + black leather). Stronger and Luxurious. (See components of the system.)
4. Steel house of main machine. Beautiful looking. (Two Choices: Silver and Black) (See components of the system.)
5. 2008 latest technologies in FIR belt. More effective and powerful. More energy -----can go up to 90°C within 5 minutes, and give out more rays

than normal belt; more safe -----100% covered with silver insulation material; belt is longer than others----- suitable for big persons. (See components of the system.)


6. Unlimited Voltage: 110V240V. Our machine can be used in any country. Transformer is in side and voltage is automatic switched. (Other machines need a switch for 110V or 220V).


7. 10 coils inside array (black and white) ----last much more times. Many companies sell 6 coils array.

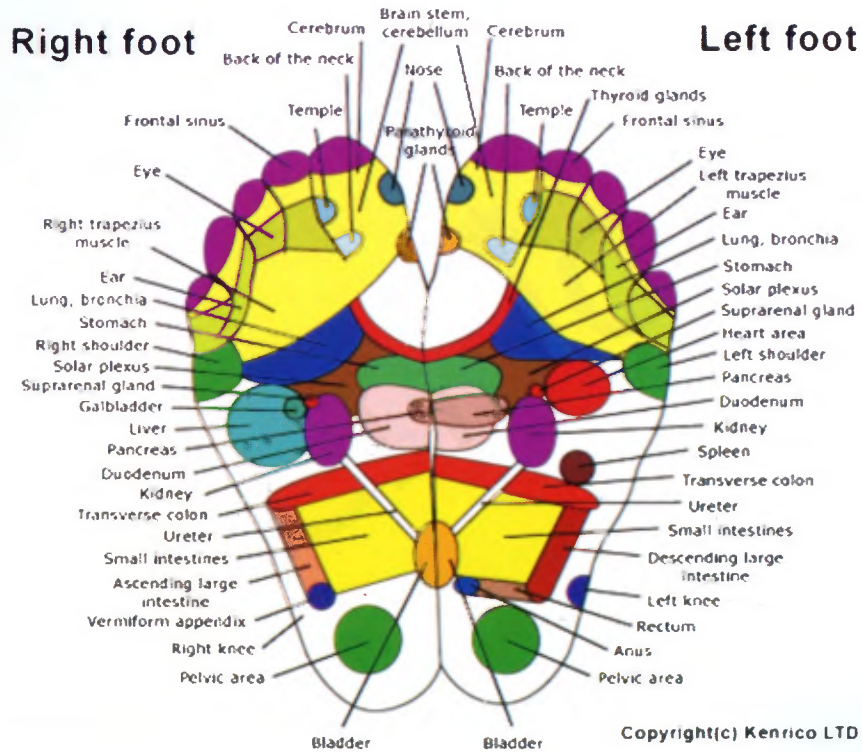

• Instruction:


1. to begin, Place a plastic bag or liner in the tub and fill with warm water, to almost completely submerge the array. Plug in the power cord, arrays, and the wristbands into the Machine. See below for control panel:

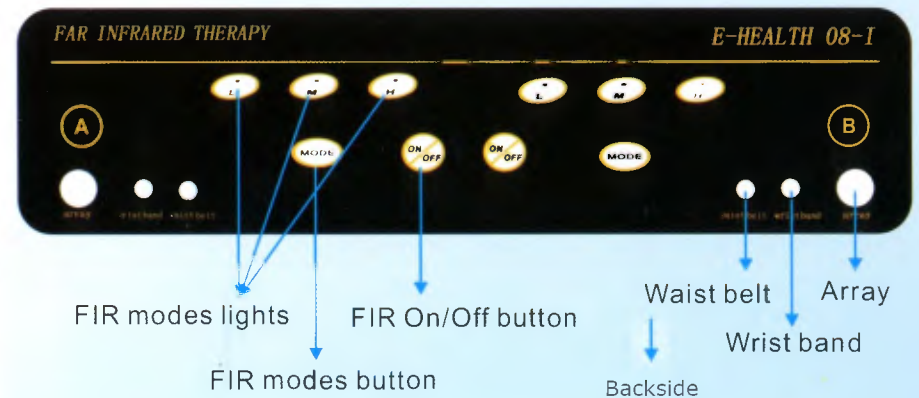

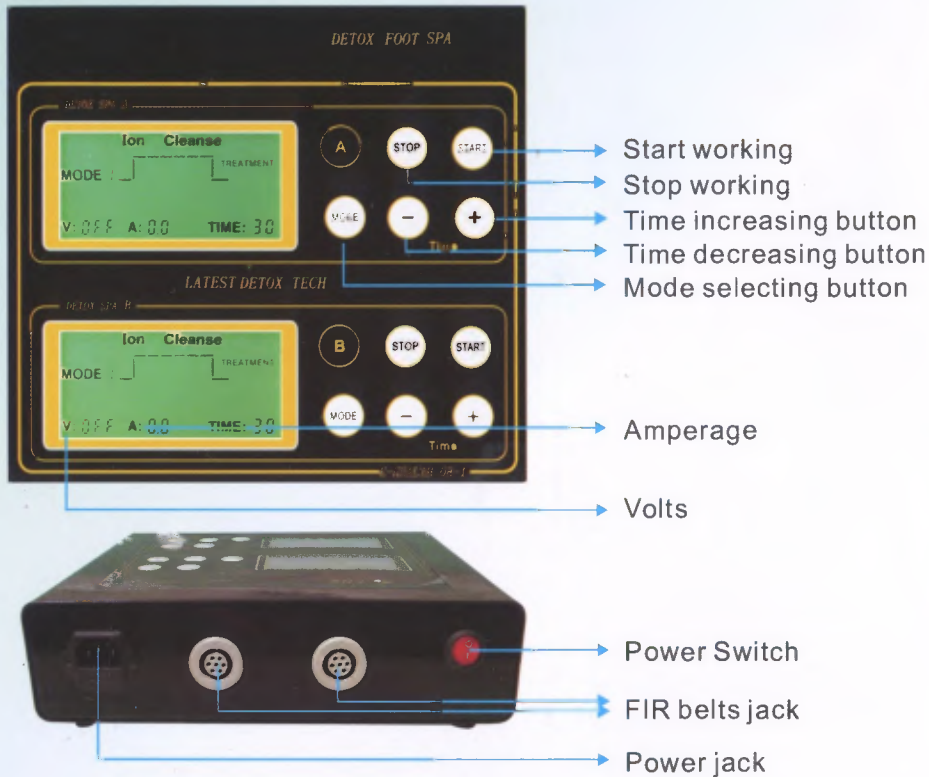

Put the wrist bands on your wrist and feet are placed in the water with the array plugged into the machine. Finally, turn on the power.

2. When you turn on the power switch, the screen shows as follows:

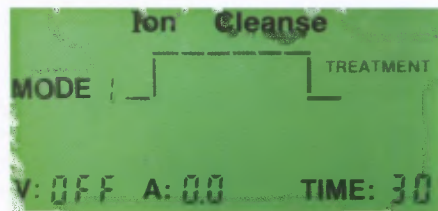

Default time is Time = 30minutes. Maximum time is 60minutes.

To adjust for more or less time, use the (+) or (-) key. 30minutes is normal set for healthy person.

3. Mode 1 is default mode. Before start therapy, Press “Mode” to select a suitable mode .There is five modes for choices .The third mode is suitable for health person.

- [1]. All positive polarity, expels acidity substance and raises PH levels.
- [2]. All negative polarity, expels heavy metals and lowers PH levels.

- [3]. (70%+)(30%-)Expels acidity and heavy metals, increases PH levels. It is the most popular, suitable for most people.
- [4]. (15+)(10-)(5+) Expels acidity and heavy metals, raises PH levels, suitable for all kinds of physicals.
- [5]. (10-)(15+)(-5) Expels heavy metals and acidity, maintain PH level, appropriate for all kinds of physicals.

4.If without salt in water, the amps show 0.6-0.9A. But Amperage should be between 1.0-2.2 for cleanse to be most effective. So you need add some salt in water. The alarm will indicate if the salt content is too high (over 2.8amps). In this instance, water must be replaced and less salt added.

5. Press “start” to begin therapy .The screen show as follow (without salt):



The machine automatically turns off when the timer ends.

Note: All water, with the exception of distilled water, has some degree of conductivity, or ability to transmit electricity. Conductivity impacts the level of ionization, which determines the effectiveness of the footbath session. At the beginning of a session, if the Footbath is reading less than 1.0 amps, add a small amount of sea salt to the water to raise the water’s conductivity. The optimum operating range is 1.0-2.2 amps.

Important note: Add as little sea salt to the water as possible, to maximize the life of the array.

• What to Expect after an ionic detox footbath Session

Common reactions from a detoxification session may include: thirst, light-headedness, mild headache, hunger, or the need to rest. Loose bowel movements may also be experienced for two to three days. Some individuals experience none of the above. Healthy individuals may feel lighter and experience a greater feeling of well being.

It is IMPORTANT that you are eating a balanced, nutritional diet, drinking plenty of water, and taking the appropriate supplements to support your unique physical/psychological needs. Remember detoxification accelerates changes in the body and healthy lifestyle choices will support this process.

The advanced Far Infrared Ray with Bamboo Charcoal Waist Belt - Instructions -

The latest studies have shown that mature bamboo charcoal (made from 5 year old bamboo from remote mountains), refined through high temperature, has a great absorption effect in purifying both air and water and adjusting humidity levels. The bamboo charcoal also emits far infrared ray, negative ions, and absorbs electromagnetic waves. The waist belt made of refined bamboo charcoal will warm your body, relieve muscle aching, increase cell activity, and increase blood circulation.

The most light we see every day is actually a combination of many different types. All lights are electro-magnetic waves. The wavelengths of lights determine their properties. Visible or invisible light can also present different colors. The Sunlight spectrum ranges from ultraviolet, on the short wavelength side, to red and infrared on the long wavelength side.

A. The eight magical efficiency of the charcoal of bamboo

1. It makes the microorganism nonpoisonous and decomposes to absorb the harmful substance
2. Supply with abundant natural mineral
3. Enable the infrared ray far warm-heat
4. Supply of the positive energy
5. Release anion prevent wearing out and purifying air
6. Covering harmful electromagnetic wave
7. It makes water soft to narrow the water Molecule
8. Deodorize and adjust the wet result

B. What is the far infrared ray

1. Biological light that can promote health

According to the life scientific research of developing rapidly in recent years, find the wavelength lies between far infrared ray and biological bearing close relation of 4- 1 microns . Form the cells of the biological organism, the principal ingredients are water and macromolecule, so it is all certain attractive to make up among molecule and molecule, atom and atom of cells, in order to reach mutual equilibrium combine, and like between the balls with the double hit of spring. And it is (such as flexible, corner, twists) to move each other with a certain fixed frequency, show the vigor of the biological life.

When the movement frequency among molecule , atom that frequency

and cell of the far infrared ray form is identical , energy will be absorbed by organism cells immediately , cause and shake the sympathetic response together , increase the vibrations in the molecule, activate and organize cells, accelerate the enzyme vigor and supply of the nutrient , promote the metabolism and health.

2. The characteristic of the far infrared ray and healthy relation

- 1) Promote the reduction reaction of the biochemistry: Reduce and corner the surplus fat, improve the physique.
- 2) Keep warm: Save the consumption of energy, keep the nutrition state of the body.
- 3) Promote the perspiration and get rid of the lower fat of skin: Reach sports, improve looks with the healthy result, but must pay attention to supplementing the moisture.
- 4) Activate physiological function and function: Activate hormone and enzyme, remove internal heavy metal, toxin, etc. and precipitate the wastes.
- 5) Invigorate: Promote the metabolism, sufficient use of nutrient component and dispel fatigue.
- 6) Physiotherapy muscles and nerve disease: Produce the positive function to the creative, lactic acid and liver; dispel the tense and not bright stuffy pain.
- 7) Inhibit the cancer cell from spreading: Reduce the destruction strength of the bad protein.
- 8) Stop the pain and diminish inflammation: Activate the immunological function.

3. Application on medical science

At present, the application in medical treatment of far infrared ray, is based on collateral channels theory of traditional medical science. And is thought to have something in common with the qigong So China's Mainland has already raised the clinical trial of the infrared ray far, the ones that were participated in are mostly Chinese medicine doctors, in physical treatment , ache application it only as for Western medicine , blood circulation not good , edema,, urine sour and higher .

It is more than one inch in depth that the far infrared ray energy can permeate in the body, its energy outputted will be transformed; make it similar to radiation energy of the body itself. So health organization can absorb, close to 93% with infrared ray that skin keep in touch, while seeing through the sleep, the far infrared ray is on the amazing function that the point gives play to of bottom of foot, promote the blood circulation of the body, this has used the concept of reflecting the treatment, so can relieve the pain muscles that the blood circulation causes.

4. The human body was made up by 70% of the moisture Molecule, and

the moisture Molecule that the plant includes is up to 90- 95% , The far infrared ray can make the moisture Molecule, often keep the resonance , sympathetic response state, produced fierce striking each other among the moisture Molecule can make volume big moisture viscount change into the very small moisture Molecule , can promote the cell activity , strengthen metabolism , achieve the goal that the human body prevents disease , even can treat persistent ailment not healed for many years . Promote the plant growing.

The far infrared ray electromagnetic wave has certain results to human health care; the following are impacts on most systems of the human body:

- 1) Blood circulatory system: If anemia, it is the high /hypotension, last apoplexy, blood obstacles competent, blood circulations angry is the bad.
- 2) Immune system: Such as pathological change of the lymphocytic system, immune system being routed.
- 3) Nerve system: Such as contain nerve lack proper care , neuralgia , have a headache , shoulder pain , pain in the back , have a backache by oneself, insomnia , physiological pain nervous all sorts of caused ache.

C. Bamboo carbon far infrared ray waistband characteristics

1. Causes the waistband to give off heat through the instrument in the DC voltage output
2. On the board is equipped with three button (low ,middle ,high) to control the waistband the temperature
3. Uses the release far infrared ray the activated carbon textile fiber for the heating element;
4. The lining has used the far infrared ray spinning and weaving product;
5. The waistband uses the high carbonization the bamboo carbon
6. Not only may warm up waist in the winter, but also has the remarkable health care function.

D. Bamboo carbon far infrared waist belts have health care function

1. The improvement waist micro cycle, invigorates the blood to remove extraverted blood, elimination lumbago, backache;
2. To treats the waist category rheumatism to have the remarkable effect;
3. To the anxious chronic nephritis, the chronic gastritis, the anxious chronic constipation and so on have the assistance treatment function;
4. Effective adjustment reproductive system, improvement month menorrhagia, the menstruation are out of balance, the premenstrual syndrome and so on;

5. The promotion blood circulation, causes fat piling up place that waist or leg department gives off heat, then the decomposition fat, reduces the size body, Healthy build.

E. Which people need the bamboo carbon far infrared waistband?

1. The acute waist sprains, the chronic strain of lumbar muscles.
2. Rheumatism cause the lumbago.
3. The anxious chronic nephritis, gastritis, the constipation, the menstruation are out of balance and so on sickness.
4. Needs to lose weight woman which likes to look good, the gentleman which the abdomen fat piles up.

Trouble Shooting

(a) If you encounter the problem of "NOT" getting any electrical current readings (THE CURRENT READS ZERO for **volts** and **amps**), the first thing confirm is DID YOU ADD SALT... as the current will **NOT REGISTER** unless SALT is added to the tub water. If you have added Sea Salt and the electrical readings are still showing ZERO then next you should confirm that the WRIST CUFF (wrist band) is on your wrist correctly allowing the Wrist Cuff metal tabs to be in contact with your skin. The two current readings are working best when they show around **1.5amps** with **14-15 volts**...

(b) Also, **MOISTEN** the Wrist Cuff CONTACTS that TOUCH to your wrist skin to elevate the Wrist Cuff's ability to read the electrical current from your wrist skin. NOTE: Similarly, I have a treadmill with a heart rate band (monitor) that goes around my chest, that if I fail to moisten the two metal tabs that touch my chest skin then sure enough, it fails to read my heart beats and shows zero on my treadmill LCD. So remember that Water Moisture helps the conductive properties of electricity.

(c) Is the WRIST CUFF plugged in the wrong place? The Wrist Cuff does "NOT" plug into the BACK of the machine it plugs in the FRONT, right beside the Array plug. The Wrist Cuff COMPLETES THE ELECTRICAL CIRCUIT (the circle of current from your body at the feet through the machine through the wrist cuff and back to the feet and again up through the machine around in a circle)... Therefore if any part of the Circuit is loose or not making PROPER CONTACT, then the LCD will show ZERO...also, no salt will cause a zero AMP reading.

(d) Please apply between a HALF and one full TABLESPOON of SEA SALT to start with, and you will notice as you sprinkle in more salt, the AMPS will rise, so add salt until you come to approximately 1.5 AMPS as the reading on the console. We have had customers tell us that ONLY Sea Salt works (you can get it at any Wal-Mart or Grocery Store in the salt & pepper

section) to further assist the conductive electricity as well as TOXIN release from your body. Sea Salt does NOT contain Iodine added but regular table salt does, so Sea Salt is best.

(e) Just to further clarify, the FIR Heating Belt (Pad) that you wrap around your stomach is what plugs into the BACK of the machine. This FIR BELT-PAD is supposed to be worn to WARM THE BLOOD in your body to improve blood flow and to further assist in releasing the TOXINS from your BLOOD and out through your feet. The FIR Heating Belt plugs into the BACK SIDE of the main Console Unit, not the front.

(f) If you happen to push the WRONG PLUG JACK into the wrong plug-in outlet hole in the Machine Console and damage the Plug Outlet Jack hole, then YOU alone are responsible, NOT US. We do not carry outlet parts, so the machine will have to be replaced if you damage it in such a way. The PLUG-IN OUTLET holes in the main Console Unit are sensitive to abuse, so if you JAM or BEND the Array Plug Jack head or Wrist Cuff plug head "TOO HARD" into their respective outlet hole you can damage the Plug Outlet Hole making the Outlet Hole loose, or damage the sensitive contacts or wires inside the Outlet Hole itself which we are NOT responsible to cover under the one year warranty. So please be gentle with the plugs and the outlet they slide into.

Warning

1. DO NOT touch power device while your HANDS ARE WET, please pay attention to your personal safety.

2. NEVER SUBMERGE the power device under water, only the Array goes in the water.

3. Children under the age of 8 **SHOULD NOT USE.**

4. Avoid striking the machine or dropping it in any way.

5. People who have low blood sugar **SHOULD EAT BEFORE USE.**

6. People who have battery operated implants (e.g. pace makers) **SHOULD NOT USE.**

7. For cleaning this device, use a clean cloth to dust the machine's surface and do not use any water or chemicals.

8. Pregnant women or nursing woman **SHOULD NEVER USE.**

9. Anyone who has received an organ transplant **SHOULD NEVER USE.**

10. People who suffer from severe heart disease that are on medication all the time **SHOULD NOT USE.**

11. People who suffer from high blood pressure that are on medication all the time **SHOULD NOT USE.**

Warranty Information

This warranty covers defects in material and workmanship for one year from the original date of delivery. We will replace any part which in our opinion is defective. If the housing has not been subjected to tampering, alteration, or improper use after delivery and has not been repaired by anyone except us. Our obligation does not include the cost of transportation. We are not responsible for damage in transit, and claims for such damage should be presented to the carrier by the customer. The housing is not warranted against misuse or use in abnormal temperature conditions.

For servicing under this warranty, return any defective part to your DISTRIBUTOR within 1 year (since the day you purchase the machine).

Warranty Registration Card

To register your warranty for this product, complete the information below.

Name

Address

City/State/Zip Code

Telephone

Date Purchased

Serial number (Base of unit)

Återförsäljare

Zon1.se

0736709495